
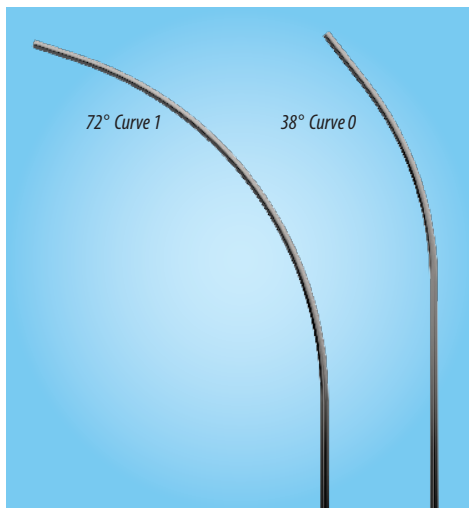


Transseptal Access without a Transseptal Needle

Specifically designed to provide support when using the SSRF™ to cross the interatrial septum

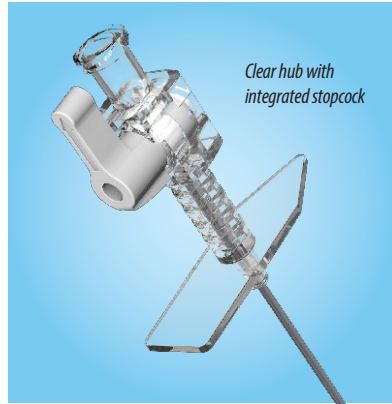

PRODUCT DESCRIPTION

The Needle Free™ Transseptal Cannula replaces a transseptal needle when using the SafeSept RF™ Transseptal Guidewire for creating the primary puncture in the interatrial septum, providing access through the septum from the right side of the heart to the left side.

Designed to be used with the SafeSept RF™, the Needle Free™ Transseptal Cannula comes in equivalent lengths of 71 cm, 89 cm, and 98 cm, and is available in Curve 0 and Curve 1 shapes.

Turn the stress that accompanies every transseptal stick into a routine, “No Hassles” procedure. The Needle Free™ Transseptal Cannula is part of a complete line of vascular access devices from Pressure Products.

FEATURES

- Clear hub to visualize devices and fluids
- Light, ergonomic hub with arrow indicator for curve direction
- Integrated stopcock with arrow indicator on handle
- Cannula provides column strength to any adult transseptal introducer system
- Smooth bore on inside of cannula to facilitate passage of the SafeSept RF™ Transseptal Guidewire
- Non-skiving design
- Mimics shape of standard transseptal needle curves
- Available in Curve 0 and Curve 1 shapes

Ordering Information

Needle Free™ Transseptal Cannula

Model	Description	Size	Length	Curve	Curve Angle
TSC071	71cm Transseptal Cannula	18G	67.7cm	0	38
TSC089	89cm Transseptal Cannula	18G	87.7cm	0	38
TSC098	98cm Transseptal Cannula	18G	96.7cm	0	38
TSC171	71cm Transseptal Cannula	18G	67.7cm	1	72
TSC189	89cm Transseptal Cannula	18G	87.7cm	1	72
TSC198	98cm Transseptal Cannula	18G	96.7cm	1	72

Packaged in quantities of 5 units per box. Individual units are not sold separately.